

## UAB "LT Biotech"

Į.K. 302303586, PVM kodas LT 100004741118,  
Rugių 21-24, LT-08419, Vilnius Reg. Nr. 127918,  
V.Į. Registrų centras Vilniaus filialas  
Tel/fax +370 5 216 02 27

### Accutase Cell Detachment Solution LTDe01

#### General Information

Accutase is a cell detachment solution of proteolytic and collagenolytic enzymes. It has been developed to meet the most demanding requirements for gentle and effective detachment of adherent cells. It is useful for the routine detachment of cells from standard tissue culture plasticware and adhesion coated plasticware.

The advantages of Accutase over the traditional Trypsin/EDTA treatment are that it is less damaging to cells leading to increased viability and no inactivation of the enzymes by inhibitors or serum is required. The Accutase solution does not contain mammalian or bacterial-derived products.


#### Applications

- Gentle detachment of a vast variety of cells including very sensitive cell lines e.g. hESCs
- Routine cell passage
- Analysis of cell surface markers
- Quiescence assays by serum starvation
- Transformation assays by oncogene transfection
- Neural crest cell migration assays, tumor cell migration assays
- Cell proliferation, apoptosis
- Production scale-up (bioreactor)
- Flow cytometry

|                        |   |
|------------------------|---|
| Appearance             | Clear, pale red to orange solution  |
| Storage and shelf life | Store at $\leq -15^{\circ}\text{C}$ Once defrosted, store at $+2$ to $+8^{\circ}\text{C}$ . Note: Do not thaw the product in a $37^{\circ}\text{C}$ water bath. After thawing, the Accutase solution may be stored for up to 2 months at $+2$ to $+8^{\circ}\text{C}$ . Do not store at room temperature. |
| Shipping conditions    | Frozen (Dry ice)  |
| Working concentration  | 1x (ready-to-use)   |
| Origin                 | Crustacea   |

#### Formulation

| Components                | mg/l |
|---------------------------|------|
| EDTA 4Na                  | 220  |
| KCl                       | 200  |
| $\text{KH}_2\text{PO}_4$  | 200  |
| NaCl                      | 8000 |
| $\text{Na}_2\text{HPO}_4$ | 1150 |
| Phenol red                | 3    |
| Accutase enzymes          | 1 x  |

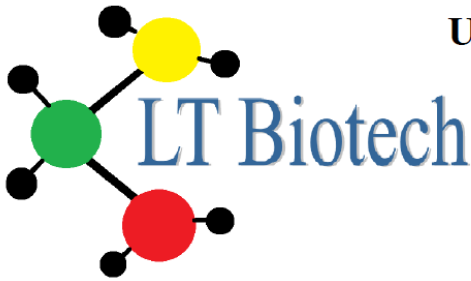

## UAB "LT Biotech"

Į.K. 302303586, PVM kodas LT 100004741118,  
Rugių 21-24, LT-08419, Vilnius Reg. Nr. 127918,  
V.Į. Registrų centras Vilniaus filialas  
Tel/fax +370 5 216 02 27

### Instruction for use

The Accutase solution is supplied as a sterile, ready-to-use, frozen liquid. Use aseptic techniques when handling this product.

1. Frozen products can either be thawed at room temperature or overnight at +2 to +8°C.  
Note: Do not thaw the product in a 37°C water bath).
2. Remove medium from culture vessel by aspiration and wash the cells with a Ca<sup>2+</sup> and Mg<sup>2+</sup>-free salt solution (e.g. Dulbecco's PBS w/o Ca & Mg; Cat. No. PBS-1A). Remove salt solution by aspiration.
3. Add Accutase to the culture dish or flask using aseptic procedures, using 10 ml for each 75 cm<sup>2</sup> of surface area. Note: Do not pre-warm Accutase in a 37°C water bath before adding to cells.
4. Return culture to 37°C incubator and allow cells to detach 5 – 10 minutes. The precise time needs to be determined for each cell line.
5. Count cells and passage as usual, no additional washes or enzyme inhibitors are required.

### Precautions and Disclaimer

This product is for research use only. Please consult the Material Safety Data Sheet for information regarding hazards and safe handling practices.